


香港地質學會

The Geological Society of Hong Kong

P.O. Box 92341, Tsim Sha Tsui Post Office, Kowloon, Hong Kong

Email: geolsoc.hongkong@gmail.com

Website: <http://www.geolsoc.org.hk>

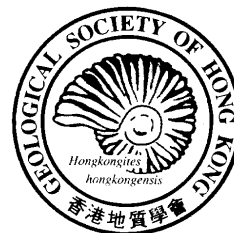

Laboratory Visit : HKUST Geotechnical Centrifuge Facility, CWB, Kowloon, Hong Kong

Date: November 23, 2013 (Saturday)

Meeting Time and Place: 9:45 am, Bus Terminus South Gate (near Shaw Studio), HKUST, Clear Water Bay, Kowloon

Trip Leader: Felix LEE

The Centrifuge is a unique and powerful tool for physical modeling of a wide range of engineering correlated problems. The facility is equipped with a 4-axis robotic manipulator, the world's first biaxial (2-D) shaking table as well as a state-of-the-art distributed data acquisition and control system. The robot can simulate construction activities such as tunneling, excavation and piling; test designs; and provide technical solutions for various geotechnical problems pertinent to Hong Kong, including short and inclined piles, loose-fill slopes and retaining walls. The shaking table is capable of simulating seismic ground response for a better understanding of earthquake-induced soil interactions such as liquefaction and deformations. Such modeling is important to regions of active seismic movement, including the Chinese Mainland. Results of the testing can be analyzed and displayed in real time by the data acquisition and control system.

Fees: Free of charge.

Transportation: Participants can reach the assembly point either from Hang Hau MTR Station or Diamond Hill MTR station or Choi Hung MTR Station by (Bus Rt 91, 91M) GMB No. 11. Or from Ngau Tak Kok MTR Station taking min greenbus 104.

[For the details of the transportation, please visit the HKUST transportation webpage . .](#)

Deadline: November 15, 2013 Friday noon

No. of participants: Due to the lab. environment and space, max. 15 quotas. Priority will be given to geotechnical field members, then members.

CPD cert.: 2 CPD hours

Enquiry: geolsoc.hongkong@gmail.com

Booking Slip (Confirmation will be emailed)

I wish to book _____ place(s) for the field trip.

Name: _____ Tel.:(office) _____ (home) _____

Mobile: _____

Email: _____

Other participants:

Name:

Member/Non-member

I and on behalf of the party listed above, accept that the Geological Society of Hong Kong and the organizers of this excursion are not under any liability for the loss of property or injury to me or any member of the party.

Signature

Date

HKUST Geotechnical Centrifuge Facility

香港科技大學岩土工程離心機實驗所

Geotechnical Centrifuge is a powerful and advanced equipment for physical modeling of various engineering related problems, such as rain-induced landslides, consolidation settlement of reclaimed lands, seismic ground response, earthquake-induced liquefaction and deformations, offshore foundation response to wave loading, seabed pipeline instability, pollutant transport in porous media, tunneling, deep excavations and many other soil-structure interaction problems under both static and dynamic loads.

岩土工程離心機是用于機體轉動時產生的離心力去加強模擬物的實際重力，這樣，很多現場的試驗都可在室內進行了。建在科大的離心機可用于模擬多种工程問題，例如暴雨產生的滑坡，填海地盤的固結沉降，地面的地震反應，地震液化及地基變形，污染物在孔隙介質中的傳播，以許多在動靜荷載下土与結構的相互作用問題等。



The 400 g-ton Geotechnical Centrifuge
400 g-噸岩土離心模型試驗機

Technical Data of Geotechnical Centrifuge:

- It is a 400 g-ton centrifuge with a nominal diameter of 8.5m (maximum acceleration =150g)
- First biaxial shaking table on centrifuge in the world for simulating earthquakes
- The second 4-axis robotic manipulator on centrifuge in the world for realistically simulating actual construction activities in-flight

離心機技術參數:

- 容量400 G.T., 最大加速度150g, 有效直徑8.5米.
- 裝有世界上首台用于地震模擬的雙向震動台.
- 裝有世界上第二台具備4軸向活動能力的機器可在離心機 運轉過程中真實地模擬土工結構的施工過程



In-flight 2D Shaking Table
離心機雙向振動台



In-flight 4-Axis Robotic Manipulator
離心機四軸機



香港科技大學

The Hong Kong University of Science and Technology

岩土工程離心機實驗所

Geotechnical Centrifuge Facility

Director: Professor Charles W.W. Ng (E-mail :charles.ng@ust.hk)
岩土工程離心機實驗所主任吳宏偉教授



HKUST Transport

Public Transport Guide

South Bus Station (New Bus Station)

<u>Route</u>	<u>From</u>		<u>Via</u>		<u>To</u>
91M	Po Lam	>	HKUST	>	Diamond Hill
91	Clear Water Bay	>	HKUST	>	Diamond Hill
91P*	HKUST	>	>	>	Choi Hung
11	Hang Hau	>	HKUST	>	Choi Hung
104 (circular)	Ngau Tau Kok	O	HKUST	O	Ngau Tau Kok

*Departing from South Bus Station, HKUST at 1745, 1800, 1815 & 1830 (only during fall & spring semesters or please refer to the announcement of KMB) while the last stop is at Choi Hung.

From Piazza via pavement next to Hang Seng Bank, please click [here](#)

From Atrium via Chia-Wei Woo Academic Concourse, please click [here](#)

North Bus Station (Existing Bus Station)

<u>Route</u>	<u>From</u>		<u>Via</u>		<u>To</u>
91M	Diamond Hill	>	HKUST	>	Po Lam
91	Diamond Hill	>	HKUST	>	Clear Water Bay
91P@	Diamond Hill	>	Choi Hung	>	HKUST
792M	Sai Kung	<>	HKUST	<>	Tseung Kwan O
11	Choi Hung	>	HKUST	>	Hang Hau
11M	Hang Hau	<>	HKUST	<>	Hang Hau
12	Sai Kung	<>	HKUST	<>	Po Lam

@ Departing from Diamond Hill PTI Bus Terminus at : 0800, 0810, 0820, 0830, 0845 & 0900 (only during fall & spring semesters or please refer to the announcement of KMB).

Shuttle Bus Service Guide

Service Period: 2 September 2013 - 20 December 2013; 4 February 2014 - 16 May 2014
(from Monday to Friday, except public holidays)

<u>Morning Schedule</u>		<u>Fares</u>
Diamond Hill (Sheung Yuen Street)	> HKUST (Piazza)	\$6.2
Time: 0815 - 0840 (every 5 - 8 mins)		
Tseung Kwan O (Po Yap Road)	> HKUST (Piazza)	\$4.5
Time: 0820 - 0845 (every 5 - 8 mins)		

Please click [here](#) for bus stop locations

<u>Afternoon Schedule</u>		<u>Fares</u>
HKUST (Piazza)	> Diamond Hill (Sheung Yuen Street)	\$6.2
HKUST (Piazza)	> Tseung Kwan O (Po Yap Road)	\$4.5
Time: 1800, 1820, 1840, 1900		

Free Shuttle Bus Service during Lunch Time

Piazza > East Point City
1235, 1245, 1255
East Point City > Piazza
1340, 1345, 1350

Shuttle bus users are required to show their student or staff ID cards upon boarding

Typhoon & Rainstorm Black Warning Arrangement

1. If typhoon signal number 8 is hoisted or rainstorm black warning is announced or has been in force after 7:00 a.m. in the morning, the University campus-bound shuttle bus service will not be provided;
2. In the case of no. 1 above , afternoon trip service departing from the University campus will provide as schedule if typhoon signal number 8 or above or rainstorm black warning is lowered later during the day;
3. If rainstorm black warning is announced during school day / office hours, the bus service will be subject to road condition and mutual arrangement between the operator and the University at the time;
4. When a warning, i.e. typhoon signal number 8 or above or rainstorm black warning has been in force but is lowered at or before 7 a.m. in the morning, the operator shall provide the service as scheduled;

University Road
大學道

North Entrance
北門入口

North Bus Station
北門巴士站

suggested outdoor route
to South Bus Station

Mr and Mrs Ho Ting Sik
Visitor Information Center
何廷禧伉儷訪客資訊中心

↓

↓

↓

Wong Cheuk She
Research Center for Environment
and Infrastructure
黃焯書科研中心

South Bus Station
南門巴士站

South Entrance
南門入口

Court 3

P

5

B1

Car Park
停車場

S H Ho Sports Hall
何善衡體育館

6

Bank
銀行

4

Bookstore
大學書店

3

The Hong Kong
Jockey Club Atrium
香港賽馬會大堂

10-12

13-15

B2

Escalators
自動電梯

17-18

16

Chia-Wei Woo
Academic Concourse
吳家瑋學術廊

25-26

21

A2

23

22

Escalators
自動電梯

27-28

29-30

24

The Hong Kong
Jockey Club
Enterprise Center
香港賽馬會
創新科技中心

33

J

K

31-32

Research &
Academic Building
科研與教學大樓

Courts A, B, C

CLP Substation

CLP Power Wind/
Wave Tunnel Facility
中電風洞實驗室

Geotechnical
Centrifuge Facility
岩土工程
離心機實驗室

Hourly Car Park
時租停車場

P

Lee Shau Kee
Business Building
李兆基商學大樓

Lee Shau Kee
Library
李兆基圖書館

Lee Shau Kee
Library
李兆基圖書館

Lee Shau Kee
Library
李兆基圖書館

Lee Shau Kee
Library
李兆基圖書館

Lee Shau Kee
Library
李兆基圖書館

Lee Shau Kee
Library
李兆基圖書館

SQ Tower 4
教職員宿舍四座

SQ Towers 5-7
教職員宿舍五
至七座

UG Hall III
學生宿舍三座

PG Hall II
研究生宿舍
二座

Stephen Kam
Chuen Cheong
Hall (PG Hall I)
張鑑泉樓
(研究生宿舍
一座)

Outdoor
Swimming
Pool
室外
游泳池

Chan Sui Kau
Chan Lam-Mo
Hall (UG Hall I)
陳瑞球林滿球
(學生宿舍七
座)

Bridge Link
連接橋

Amphitheater
圓形露天劇場

UG Hall I
學生宿舍
一座

Jockey Club
賽馬會樓

(UG Hall VI)
(學生宿舍六
座)

S H Ho Tower
何善衡樓

Unilodge (South Block)
水雲軒(南座)

Tin Ka Ping Hall
田家炳廳

Lo Ka Chung University Center
盧家驊大學中心

Towers C & D
C及D座

Staff Quarters (SQ) Towers
教職員宿舍八至十一座

Court 7

Lee Shau Kee
Campus
李兆基校園

To Tseung Kwan O & Clear Water Bay
往將軍澳及清水灣



THE HONG KONG GEOTECHNICAL CENTRIFUGE AND ITS UNIQUE CAPABILITIES

C.W.W. Ng, P. Van Laak, W.H. Tang, X.S. Li and C.K. Shen

Department of Civil Engineering, the Hong Kong University of Science and Technology
Clear Water Bay, Kowloon, Hong Kong

ABSTRACT

Recently there has been rapid development in geotechnical centrifuge modelling technology world-wide, and centrifuge testing is now commonly used for reduced-scale physical modelling of geotechnical and geo-environmental systems. In this paper, a new 400 g-ton state-of-the-art geotechnical centrifuge facility in Hong Kong is introduced and some of its unique and advanced capabilities such as in-flight bi-axial shaker, 4-axis robotic manipulator, and data acquisition and control system are described. This facility is dedicated to serve not only the University but also the entire geotechnical community in Hong Kong and around the world.

KEYWORDS: Geotechnical centrifuge, bi-axial shaker, robotic manipulator, geophysical modelling

1. Introduction

As illustrated in Figure 1, geotechnical centrifuge modelling has now become a widely-used modelling tool to complement more conventional numerical analysis and field testing or monitoring (Schofield, 1980 & 1998; Ko, 1988; Van Laak et al., 1994; Ng & Springman, 1994; Ellis et al., 1995; Ng et al., 1998a & b; Kimura, 1998). Each approach has its own advantages in terms of quality of result, time and cost. Particularly in cases where there are uncertainties in the applicability of a proposed design methodology, use of more than one approach permits calibration of results against each other and verification of conclusions drawn.

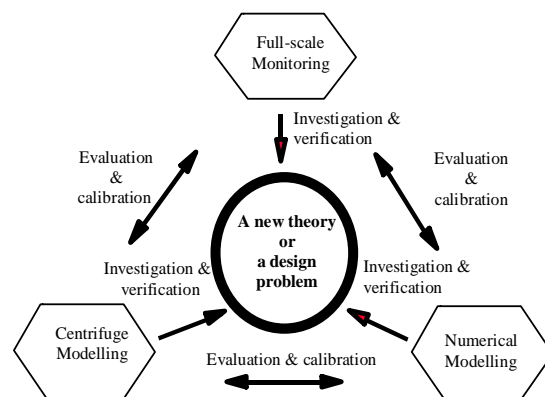


Figure 1 The inter-relationships between field monitoring, numerical and centrifuge modelling

2. Basic principle of geotechnical centrifuge modelling

For most applications in civil engineering, the major and fundamental

difference in engineering behaviour between soil and steel (or concrete) is that the behaviour of the former material is strongly dependent on stress state and stress path, whereas the latter is generally not. This implies that any laboratory or field test involving soils which fails to correctly model stress state and stress path will likely produce incorrect results. In particular, small-scale geotechnical model tests conducted under one gravity (9.81 m/s^2) can be very misleading because the stress field within the prototype is not modelled accurately, leading to incorrect dilative instead of contractive behaviour at low stresses. How to create a correct stress field effectively and economically for physical modelling of geotechnical problems has been a major challenge for engineers and researchers for many years.

The basic principle of centrifuge modelling is to recreate the stress conditions which would exist in a full-scale construction (prototype), using a model on a greatly reduced scale. This is done by subjecting the model components to an enhanced body force, which is provided by a centripetal acceleration of magnitude ng , where g is the acceleration due to the Earth gravity (i.e. 9.81 m/s^2). Stress replication in an n th scale model is achieved when the imposed "gravitational" acceleration is equal to ng . Thus, a centrifuge is suitable for modelling stress-dependent problems. Moreover, reduction of time for model tests such as consolidation time can be achieved by using a reduced size model.

Figure 2 shows a plan view of a model earth dam rotating at a constant angular velocity ($d\theta/dt$) in a centrifuge test. For a typical static model test, $r(d\theta/dt)^2$ provides the artificial "gravitational" acceleration ng , which is used to increase the body force of the model.

2.1 Scaling Laws

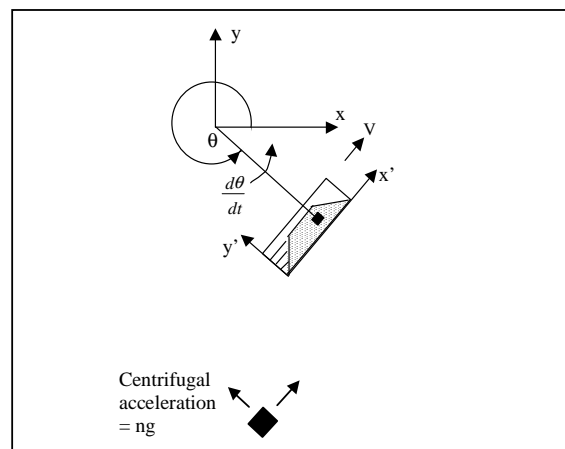


Figure 2 A plan view of model dam in a centrifuge test

For centrifuge model tests, model laws are generally derived through dimensional analysis, from the governing equations for a phenomenon, or from the principles of mechanical similarity between a model and a prototype. Some commonly used scaling laws are summarised in Table. 1.

It can be readily deduced from Table 1 that the stress level of a 100m tall slope can be correctly modelled by using a 1-metre slope model when it is subjected to an elevated "gravitational" acceleration of $100g$ (i.e., $n=100$). Also, a four-hour centrifuge modelling at $100g$ can correctly simulate a prototype soil settlement problem consolidated for more than 4.5 years (i.e., $4xn^2$ or 4×100^2 hours). Substantial time reduction and hence cost savings can be achieved by adopting the centrifuge modelling technique.

3. Development of the Hong Kong Geotechnical Centrifuge Facility

One of the largest and most advanced geotechnical centrifuges in the world is being commissioned on the campus of the Hong Kong University of Science & Technology (HKUST). Figure 3 shows the comparisons of capacity between the HKUST centrifuge and others in the world.

Table. 1 Some common scaling laws for centrifuge tests

Parameter	Scale (model/prototype)
Acceleration	n
Linear dimension	1/n
Area dimension	1/n ²
Volume dimension	1/n ³
Stress	1
Strain	1
Mass	1/n ³
Density	1
Unit weight	n
Force	1/n ²
Bending Moment	1/n ³
Bending Moment / unit width	1/n ²
Flexural stiffness / unit width	1/n ³
Time (dynamic)	1/n
Time (consolidation/ diffusion)	1/n ²
Time (creep)	1
Pore fluid velocity	n
Concentration	1
Velocity (dynamic)	1
Frequency	n




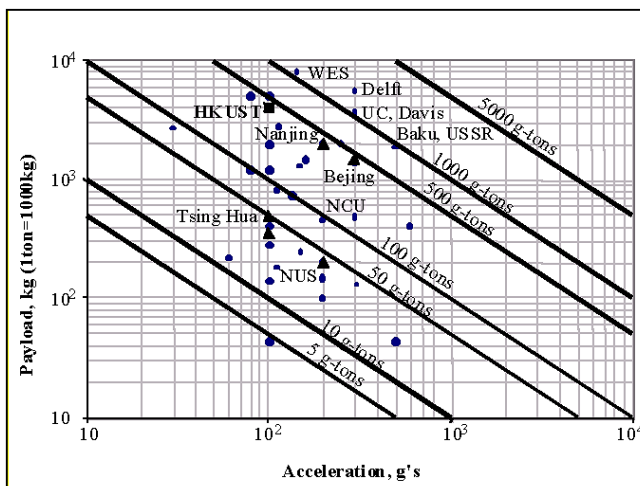


Figure 3. Capacities of major geotechnical centrifuges in the world

The centrifuge is installed in a new building, located a short distance from the main academic building on campus. The centrifuge facility has a total of 255 m² of

office and general laboratory space. In the main laboratory area, a 20 tonne capacity overhead gantry crane is available to move the pre-cast concrete panels above the centrifuge enclosure and to load and unload the centrifuge model containers. The crane is also used to interchange the static platform and shaker when required. The unused platform or shaker is stored in a recess in the floor of the centrifuge enclosure. The centrifuge is monitored using CCTV cameras and microphones, and an intercom is used to communicate between the centrifuge chamber and control room during model checkout. The hydraulic power supply is located below the main laboratory area in a room adjacent to the centrifuge.

Figure 6 shows details of an elevation view of the centrifuge. The centrifuge arm is supported on a vertical drive shaft running on a pair of pre-loaded tapered

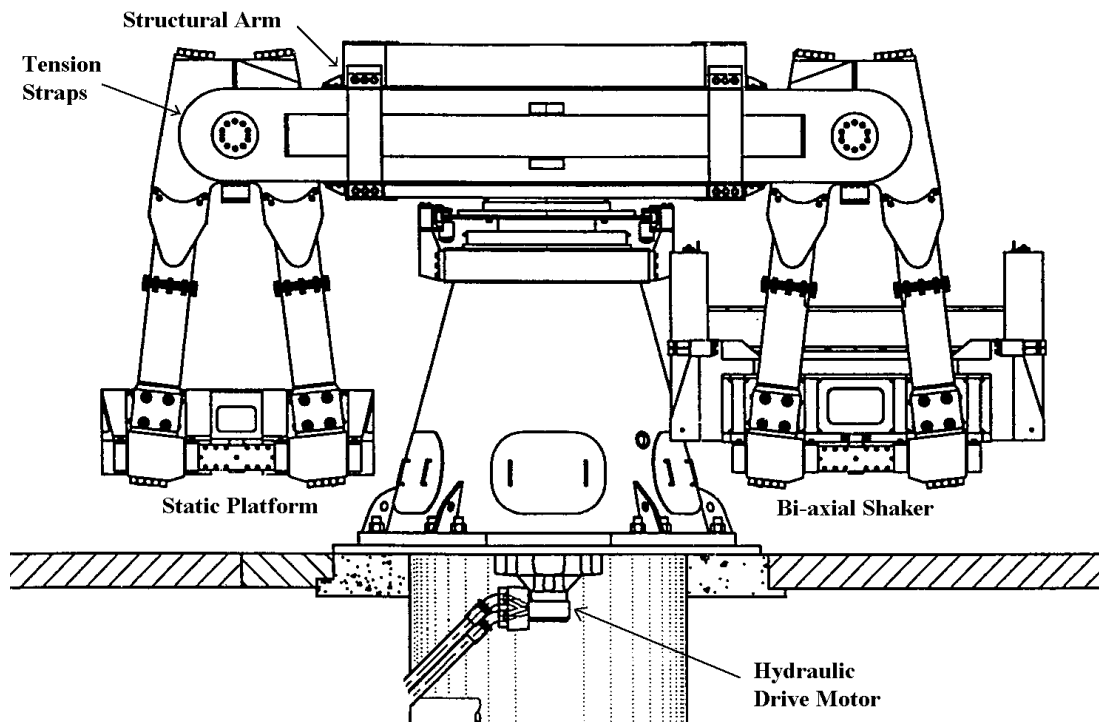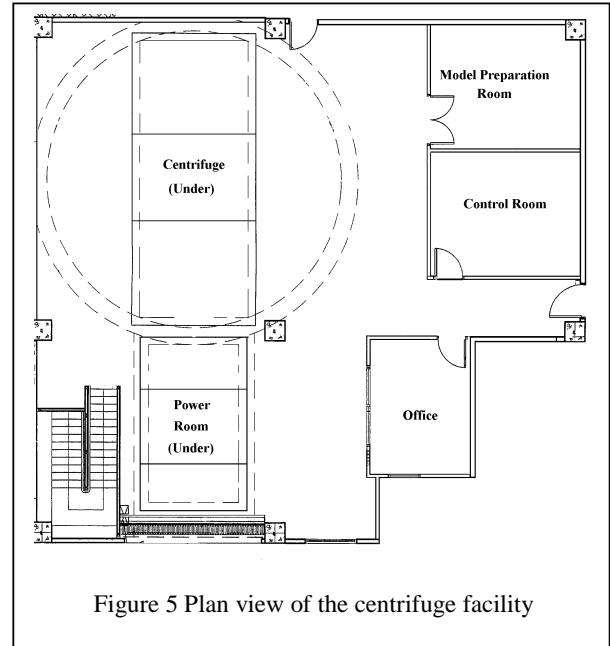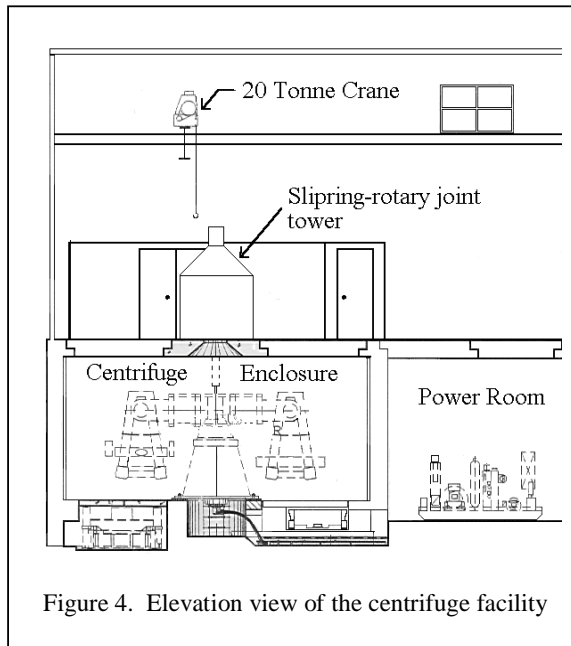

roller bearings. It is driven by a hydraulic radial piston motor directly coupled to the lower end of the vertical drive shaft. Balancing of the centrifuge is accomplished using weights of various sizes placed on the platform not used for

the model container. Careful bookkeeping is used to ensure the centrifuge is closely balanced to maintain stresses within structural limits and to prolong bearing life.

The 400 g-ton centrifuge has a rotating arm of 3.4 metres nominal radius. In total, three swinging platforms have been manufactured. Two platforms are identical and are designed for non-shaking tests. Each of these 'static platforms' can accommodate a model of up to 1.5m x 1.5m x 1m in size and up to 40,000 N in weight. The third platform comprises the bi-axial shaker and associated structural supports, hydraulic manifolds and reaction mass. The shaker slip-table can accommodate payloads of up to 0.6m x 0.6m x 0.4m and up to 3000 N in weight.

Table 2 summaries some key specifications for the centrifuge. Note that for static tests, the centrifuge can be operated at up to 150g whereas for dynamic tests, the bi-axial shaker is designed to operate at up to 75g.

The centrifuge drive uses a variable volume pressure compensated pump driven by 150 kW electric motor. A computer is used to control the centrifuge and the shaker. Besides controlling the speed of the centrifuge, the computer is used to monitor operational parameters such as imbalance forces in the arm, temperature of the main bearings, hydraulic fluid temperature and pressure, and status of safety interlocks. The computer is also used to implement the sequence of valving operations required for operating the shaker system during shaking tests. Input data and measured parameters from each run are automatically logged to the computer hard disk to facilitate long-term monitoring of machine performance and scheduling of routine maintenance.

A fluid rotary joint is provided for supplying pressurized air, water, and hydraulic fluid to instrumentation mounted on the centrifuge during spinning. An electrical slip ring assembly is provided for transmitting electrical power and signals to and from centrifuge instrumentation, for control of various devices and acquisition of experimental

data. Some specifications for the assemblies are given in Table 3.

4. The bi-axial shaker

The HKUST centrifuge incorporates a bi-axial servo-hydraulic shaker, to be used for simulated seismic excitation. In consideration of the facts that earthquake motions are multi-directional in nature and that many uni-directional centrifuge earthquake simulators are already available, HKUST decided to develop a bi-axial shaker in order to simulate earthquake motions in two horizontal directions simultaneously.

Because large shaking forces in two directions are possible with this shaker, development of the centrifuge and the shaker was carried out simultaneously, with the shaker designed as an integral part of the centrifuge. This integrated approach to the design was adopted in order to isolate the shaking forces from the centrifuge to as large an extent as possible and to produce high quality shaking motions.

To facilitate installation and maintenance, and to permit operation of the centrifuge at accelerations greater than 75 g for static tests, the shake table and its bucket form a single assembly that is removable from the suspension arms and replaceable by an optional static bucket.

With a payload weight of 3000N, the total moving weight (payload and shake table hardware) is about 10,000N. To optimize the dynamic behaviour of the in-flight shaker, a large reaction mass (4000 kg) has been incorporated into the design. The shake table is supported by hydrostatic self-aligning pad bearings, and utilizes two pairs of servo-actuators for each of the shaking directions (one pair in the tangential direction of the centrifuge rotation and the other in the direction of bucket swing-up). Each pair of actuators is located on opposite sides of the shaking platform and corresponding pairs are

Table 2. Technical specifications for the centrifuge.

Key item	Specification
Payload capacity	400 g-tons
Arm radius	4.2 m to the base of the swinging platform
Maximum acceleration	150g (Static tests) 75g (Dynamic tests)
Payload size	1.5mx1.5mx1m for static tests, 0.6mx0.6x0.4m for dynamic tests

Table 3. Specifications of slip ring and rotary joint assembly.

Key item	Specification
Slip rings	32 for analog signals, 8 for analog return, 16 for power
Co-axial cable channels	8 for video, and high frequency equipment, 4 high quality for digital signals (computer network)
Air ports	2 at 700 kPa, 0.05 m ³ /min
Water ports	2 at 1400 kPa, 40 liters/min

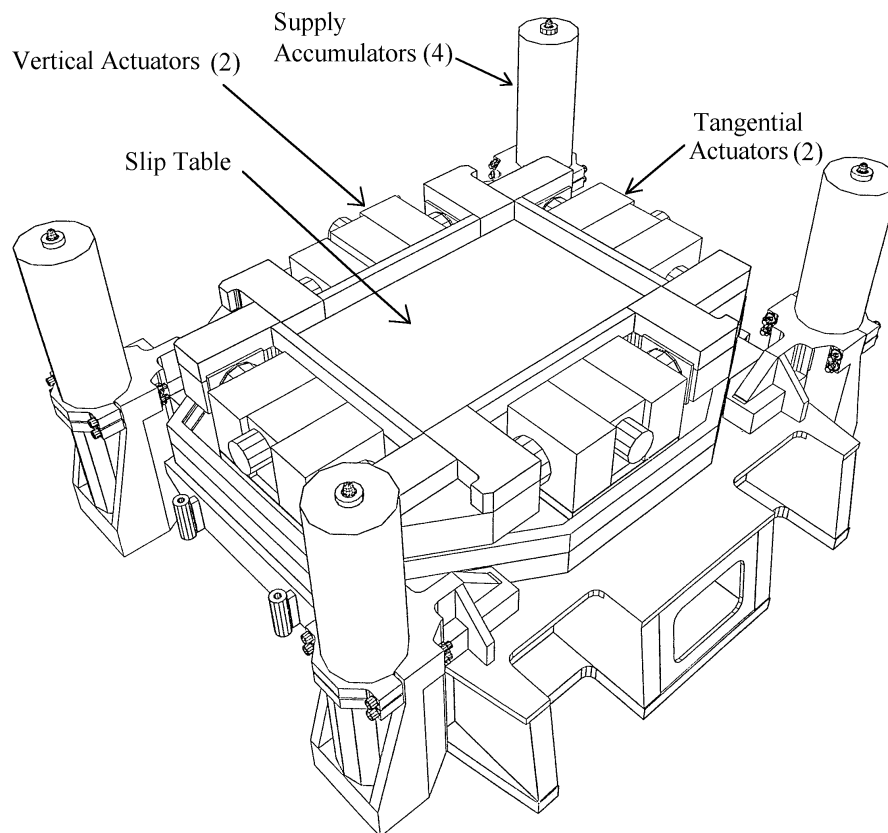


Figure 7 The bi-axial shaker

designed to act as a unit, applying identical forces to each side of the slip table. The

motion of the slip table is then ostensibly a superposition of translations in the two

orthogonal directions, with no rotations in the plane of shaking.

The shake table receives oil from a 10 gpm, 35 Mpa variable volume pressure compensated pump. This pump supplies oil to the shaker through a 35 Mpa hydraulic rotary joint mounted near the top of the centrifuge. Each pair applies a force of 32,000 kgf per axis. For the moving mass of just under 1000 kg, this force yields 35 g of acceleration in each shaking direction, and for a typical seismic signal, a shaking duration of 2 seconds for each accumulator charging is achievable. Several key technical specifications of the shaker are listed in Table 4.

To achieve the target shaking motions, a combined analog and digital control system is used. The analog portion of the controller is used primarily to eliminate the actuator redundancy by using pressure feedback to prevent actuator 'fighting'. The digital shaking control system is implemented using a computer equipped with a high speed analog/digital interface. By using an error correction algorithm, the digital controller can accurately generate a variety of input shaking motions, including simple harmonic motions and scaled earthquake motions.

5. The 4-axis robotic manipulator

An advanced 4-axis robotic manipulator (second of its kind in the world) has been developed for simulating details of construction activities such as soil nailing, pile driving, tunnelling and excavation in-flight (Figure 8). The robot incorporates a standard tool changer to permit interchanging tools without need for stopping the centrifuge. It can operate either in accordance with a sequence of pre-programmed instructions, or in response to real-time commands from an operator, in a "fly-by-wire" mode of operation. Key specifications are listed in Table 5.

6. Data Acquisition and Control System

A new state-of-the-art distributed data acquisition and control system has been developed to support sophisticated modelling activities. The system uses multiple networked data acquisition servers to provide high-speed real time data acquisition on a large number of input channels. Custom-designed computer-controlled signal conditioning provides adjustable amplification of low-level signals from measurement transducers at the source, providing signals of very high signal-to-noise ratio and exceptionally high quality. A dedicated file server is used to archive experimental data. In addition, a relational database is used to store all information required to accurately interpret data acquired during testing.

An additional benefit of using a networked system is that it is possible to make test data available to researchers in widely dispersed geographical locations over the internet via a specially designed web page. This capability, which is currently being developed at HKUST, will make it possible for clients to view experimental data in nearly real-time and to interact with the centrifuge operator while a test is in progress. Once the test is finished, the complete experimental data set can be accessed remotely via the web page for further analysis.

7. Ancillary equipment

Several model containers have been developed. One is a rectangular aluminum container, designed for static tests, having inside dimensions of 1270mm x 1245mm x 850mm, and incorporating an optional acrylic window. The container has been designed to limit average strains within the model to 0.025% during spin-up. For dynamic tests, a laminar-type flexible walled container has been developed. This container is cylindrical in shape, with an inside diameter of 550mm and a height of 500mm. The container is constructed of fifty-two aluminum rings, with roller

Table 4. Technical specifications for the shaker.

Key item	Specification
Shaking direction	Two prototype horizontal directions
Maximum shaking acceleration	35g
Maximum shaking velocity	750 mm/sec
Shaking frequency	0-350 Hz

Table 5 Key specifications for the HKUST Robotic Manipulator

Key item	X-axis	Y-axis	Z-axis	θ -axis
Stroke	1.008m	0.839m	0.305m	270 °
Maximum Speed	6.67 cm/s	6.67 cm/s	3.47 cm/s	10 °/s
Accuracy	1.0 mm	1.0 mm	1.0 mm	1.0 °
Load Capacity	1000N	1000N	5000N	5 N-m

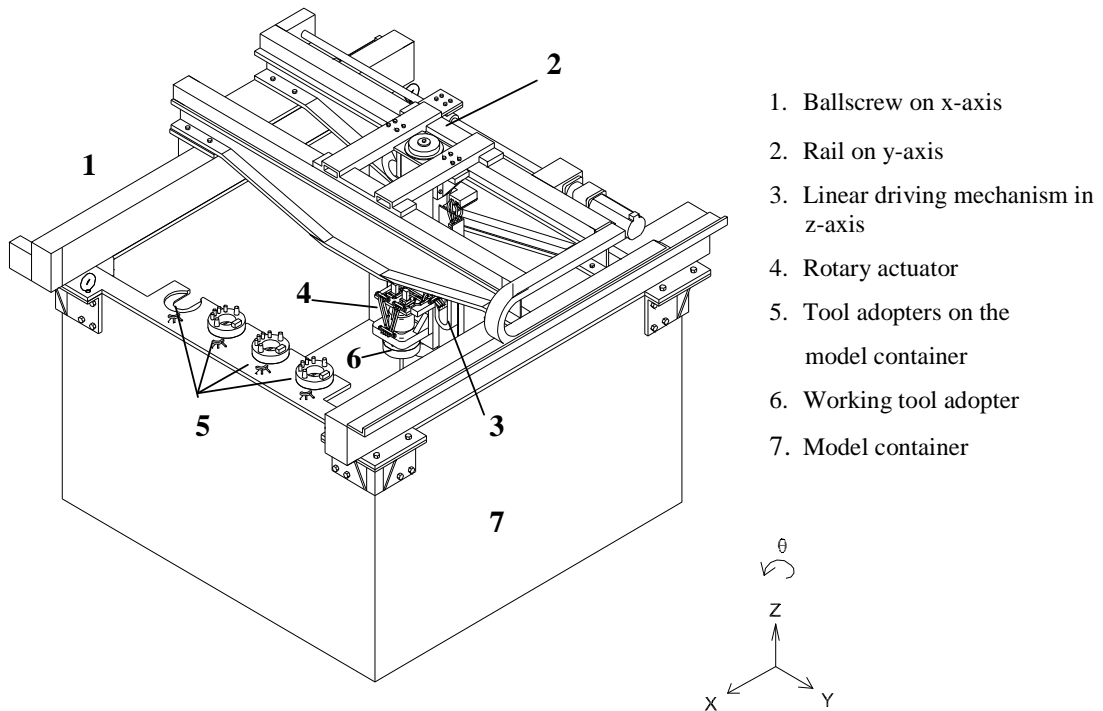


Figure 8 The 4-axis robotic manipulator

bearings used to provide the low-friction interface between adjacent rings. The laminar container has been designed to provide the necessary model boundary conditions for bi-directional shaking.

8. Summary

The state-of-the-art 400 g-ton geotechnical centrifuge at the Hong Kong University of Science and Technology is equipped with a unique bi-axial shaker, an advanced robotic 4-axis robotic manipulator and a modern data acquisition and control system. It represents a unique research facility that will be used to help academia and the construction industry advance their

knowledge and technology in the 21st century.

References

- Ellis, E.A., Ng, C.W.W. & Springman, S.M. (1995), "Piled bridge abutments on soft clay", *Proc. 10th Asian Regional Conf. on Soil Mech. & Fdn. Engng, Beijing. Vol. 1*, pp. 193-196.
- Ko, H.Y. (1988), "Summary of the state-of-the-art in centrifuge model testing", *Centrifuge in Soil Mechanics*, pp. 11-28.
- Kimura, T. (1998), "Development of Geotechnical centrifuge in Japan", *Proc. Centrifuge 98, Tokyo, Pre-print volume*, pp. 23-32.
- Ng, C.W.W. & Springman, S.M. (1994), "Uplift resistance of buried pipelines in granular materials", *Proc. Centrifuge 94, Singapore*, pp. 753-758.
- Ng, C.W.W., Springman, S.M. & Norrish, A.R.M. (1998a), "Centrifuge modelling of spread-base integral bridge abutments", *Journal of Geotechnical and Environmental Engineering, ASCE, Vol. 124, No. 5*, pp. 376-388.
- Schofield, A.N. (1980), "Cambridge Geotechnical Centrifuge Operations", *Géotechnique*, **30**, No. 3, pp. 227-268.
- Schofield, A.N. (1998), "Geotechnical centrifuge development can correct a soil mechanic error", *Proc. Centrifuge 98, Tokyo, Preprint volume*, pp. 1-8.
- Shen, C.K., Li X. S., Ng, C.W.W., Van Laak, P.A., Kutter, B.L., Cappel, K. and Tauscher, R. C. (1998), "Development of a Geotechnical Centrifuge in Hong Kong", *Proc. Centrifuge 98, Tokyo*, pp. 13-18.
- Van Laak, P.A., Elgamal, A.W. & Dobry, R. (1994), "Design and performance of an electro-hydraulic shaker for the RPI centrifuge", *Proc. Centrifuge 94, Singapore*, pp. 139-144.

Acknowledgements

The authors wish to thank the Central Allocation Vote of the University Grants Committee and the Hong Kong University of Science and Technology for providing funds to construct the geotechnical centrifuge.